




# HIFU 4D

High Intensity  
Focused Ultrasound

- Face Lifting ■
- Skin Tightening ■
- Wrinkle Removal ■
- Skin Rejuvenation ■
- Body Slimming ■
- Vaginal Tightening ■

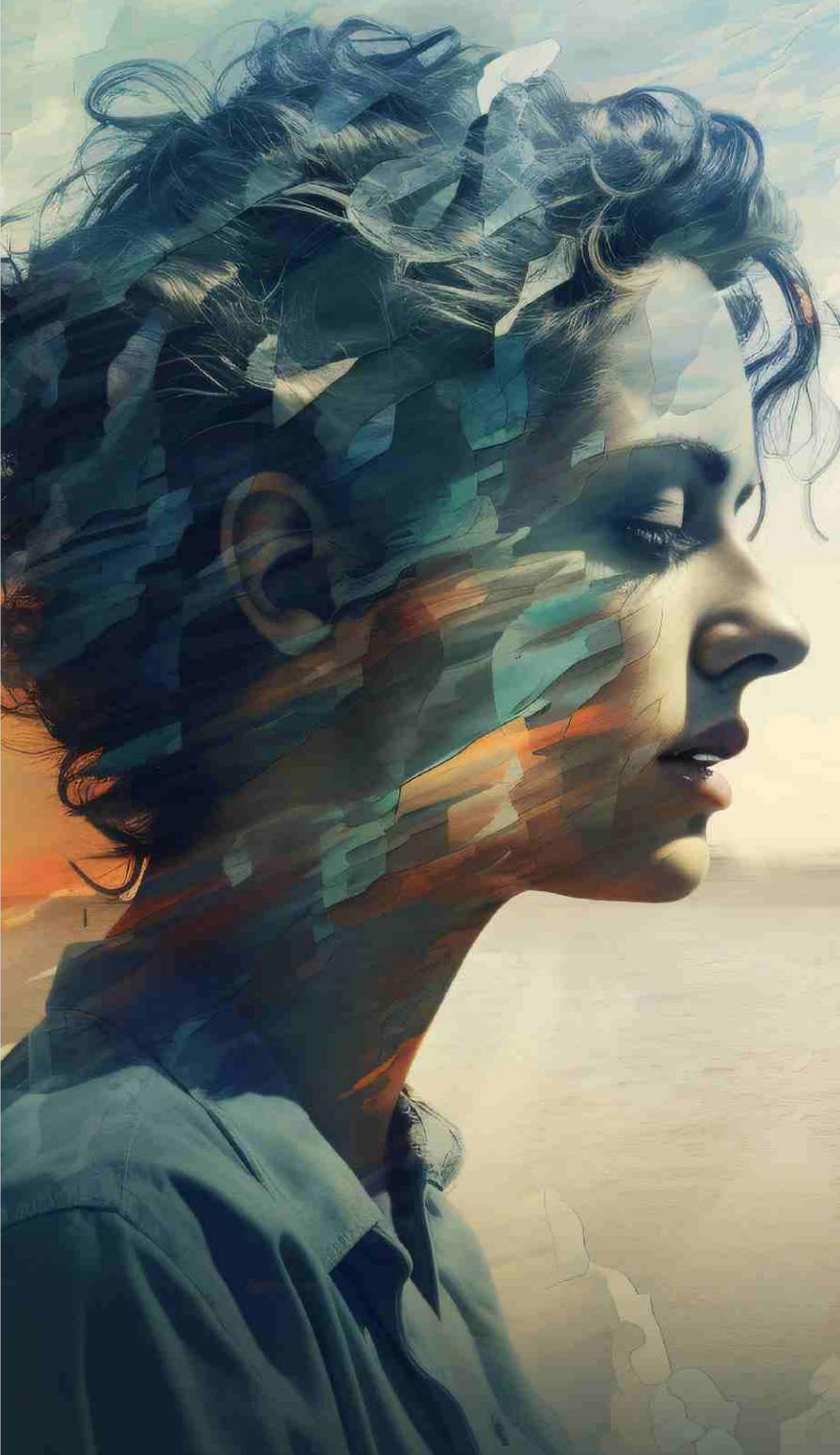

# HIFU 4D

## High Intensity Focused Ultrasound

HIFU (High-Intensity Focused Ultrasound) treatments use focused ultrasound energy to penetrate deep into the skin and target the SMAS layer, the same layer tightened during a surgical facelift. This energy generates heat, causing tissue contraction and triggering the body's natural healing process, boosting collagen and elastin production, leading to firmer, lifted skin. Here's a more detailed breakdown of how HIFU works:

**The safest & most affective HIFU for non-surgical face & body lifting and tightening with no downtime and lasting results.**




## HIFU - Handles




## HIFU - Treatment Theory


HIFU is the art-of-state technology for non-surgical skin lifting and body shaping, the full name is High Intensity Focused Ultrasound. HIFU directly delivers heat energy to skin and subcutaneous tissue that can stimulate and renew the skin's collagen and thus consequently improving the texture and reducing sagging of the skin. HIFU will heat the SMAS up to  $60^{\circ}\text{C} \sim 70^{\circ}\text{C}$ , the heat energy raises thermal coagulation point which accelerates skin rejuvenation process and maximizes tightening effect of the skin layer. Over 90-180 days, the wound-healing response stimulates long-term tissue remodeling and leads to further lifting and tightening, with results that can last 1 to 3 years.

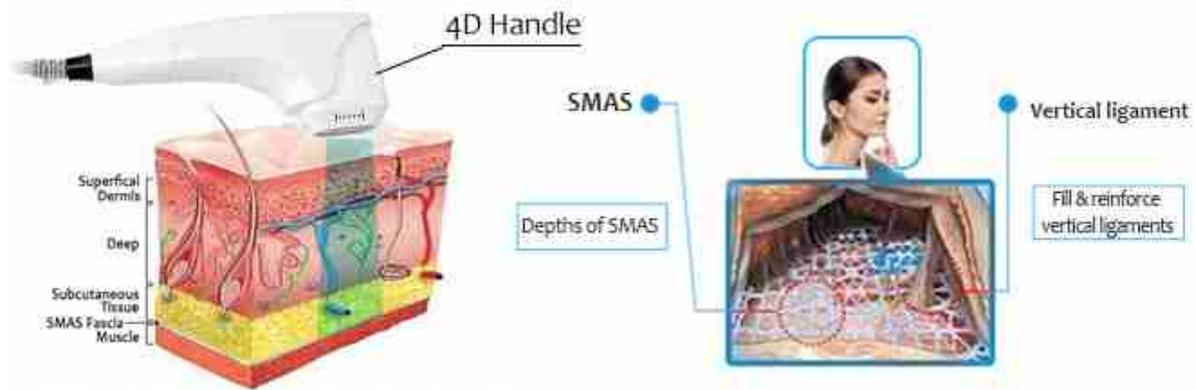

HIFU acts on vaginal mucosa and muscular tissue, generates extensive and regular thermal effect, gets the instant tightening and lifting result. 4.5mm targets on the SMAS layer, which will lift and tighten the vagina. 3.0mm targets on the deep dermis, which will improve the vagina lasting elasticity. These peeling channels will stimulate massive fibrocytes regeneration. The patented comfort technology ensures non-invasive and safety.

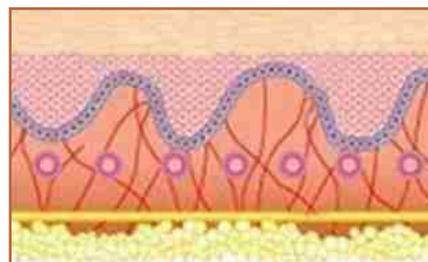

Cartridge  
D4.5  
4-7MHz

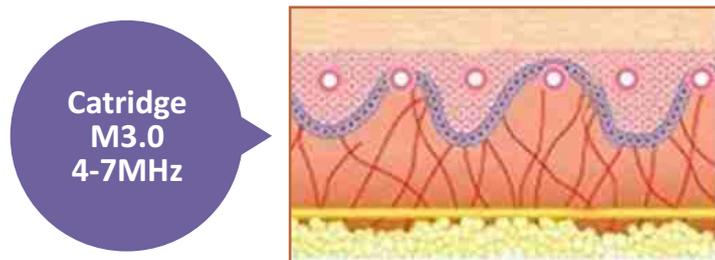

Cartridge  
M3.0  
4-7MHz

## Features & Advantages



1. Two modes, fast or slow mode can be switched.
2. Accurately affects different depths of skin.
3. Thermal stimulation on the fat & fascia layer(SMAS).
4. The operation is simple and convenient.
5. No need for consumables.
6. Immediate effect after treatment.
7. The effect can last 18-24 months at a time.
8. One shot can be output 12 lines.

## Applications

### Lifting & Tightening

Face lifting and skin tightening forehead, face and neck



### Active Collagen

Recover skin elasticity, anti-aging.





### Lift Laxity Skin

Tighten skin laxity around the eyes, cheeks, skin, chin and other parts.



### All kinds of Wrinkles

Deep wrinkles like forehead lines, crow's feet, lip and frown wrinkles.

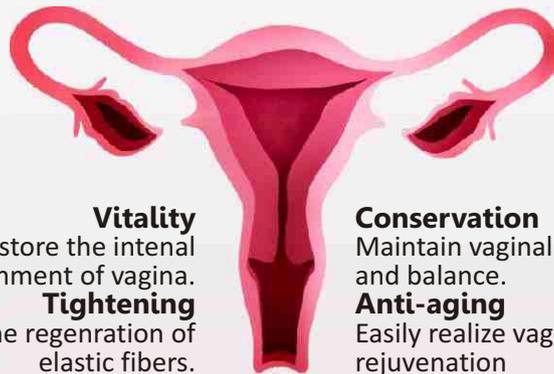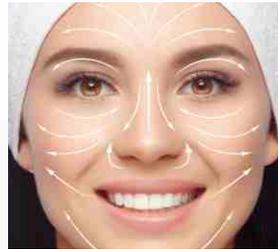

### Vitality

Restore the internal environment of vagina.

### Tightening

Promote the regeneration of elastic fibers.

### Conservation

Maintain vaginal health and balance.

### Anti-aging

Easily realize vaginal rejuvenation

## Specification

Model Name: **HIFU-V-4D** | Handles : **4D Handle, V-max Handle, Vaginal Handle (Optional)**

4D Cartridge : **3.0mm, 4.5mm (Optional:16mm, 13mm, 10mm, 8mm, 6mm)**

V-max Cartridge : **1.5mm, 3.0mm, 4.5mm (Optional: 8mm, 13mm)** | Vaginal Cartridge : **3.0mm, 4.5mm**

Output Frequency : **4-7MHz** | Output Voltage: **10-200 W** | Pack Dimension: **20 X 40 X 54 cm**

Gross Weight : **11 kgs (Without the base)** | Power Supply : **AC110/220V±10%, 50 Hz-60Hz**



+91-94460-01470  
+91-94460-01472



info@evolutionmed.com  
www.evolutionmed.com