



# 360° Cryolipolysis Cryocool



The Five Applicator  
Body Shaping System

# 360° Cryolipolysis Cryocool



## Treatment Theory


After the adipose tissue of the target care site is sucked into the cooling probe, the cooling probe releases a 360-degree surrounding cooling energy to act on the fat cells without damaging any surrounding tissue. Cryocool's steady and powerful cooling energy triggers a natural cell death process called "apoptosis", after which apoptotic fat cells are excreted through the body's natural metabolic process. Cryocool targets fat cells for permanent removal from the body, resulting in significant fat loss in the treated area.



Many of us have stubborn fat despite diet and exercise.



Warm-up technology to accelerate blood circulation.

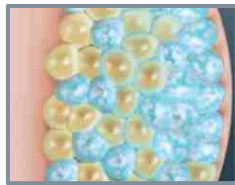

Applicator 360 degrees around fat cells without damaging any surroundingsue.



In the weeks to follow, your body naturally process the fat and eliminates these dead cells.



Cryocool procedure results are long-term, as treated fat cells are gone for good.

## Technology Advantages



Four Handles Working Together



5 applicator easy replace



360° cooling, effective in 30 mins




Medical silicone ring Safe & Comfortable



Medical filter cotton safe & convenient




Humanized small screen design, easy to operate

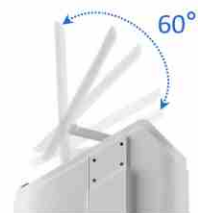


Folding screen to adapt to different scenarios

## 360° Multi Direction Multilevel Sterioscopic

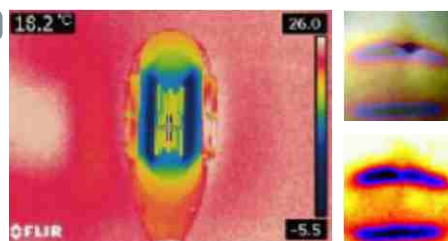

Cooling on Both Sides

Traditional Cryolipolysis by two Chilling Plates, Only 40% effective for treatment area.

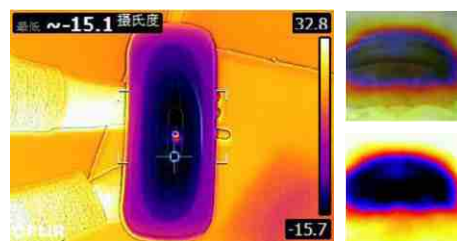


360° Refrigeration

NEW 360° cooling applicator. 100% effective for treatment area. The lowest temperature can reach -15°

**Handle  
Screen for  
Real-time  
Control.**




**Detachable  
Treatment  
Heads**

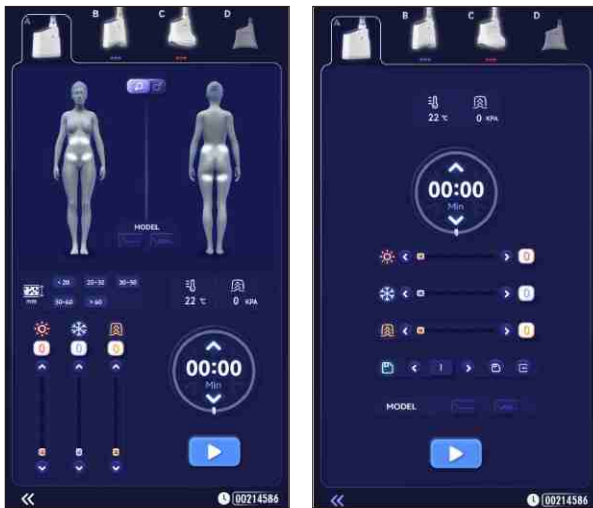

**Before / After**




**Five Applicators, More Treatment  
Areas & More Effective**

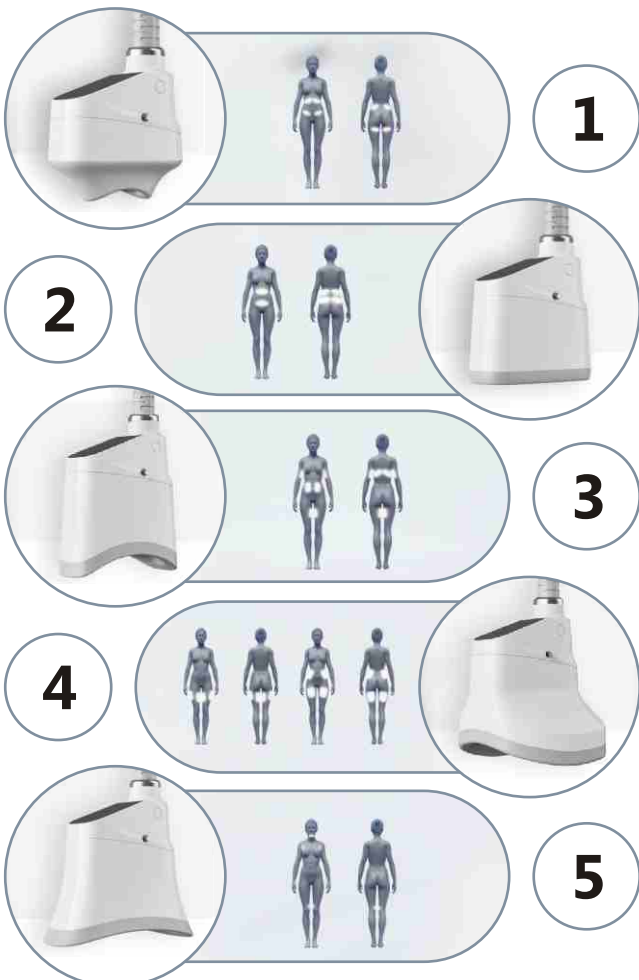


## Who is an Ideal Candidate for Cryocool?


Men and women who want to achieve cellulite reduction and inch loss without surgery or downtime and who are at a healthy body weight with no medical contraindications may be ideal candidates for 5 Applicator Cryocool.

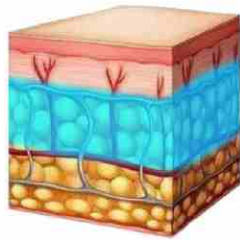

First Day

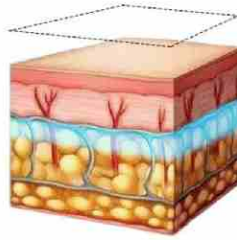

30 Days

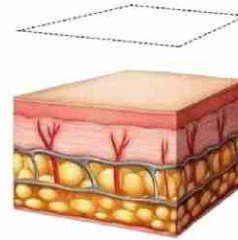

90 Days

*A course of treatment is 6 times, each time only takes 30 minutes.  
Do it at least 1 times one month and 3 months in a row easy and fast.*

## Why Chose 360° Cryolipolysis Cryocool?



Double Chin | Under Jawline | Bra Fat | Upper Arm | Abdomen/Side  
Thigh | Back Fat | Waist/Back | Flank | Banana Roll

## Specification

Input Power : **1200W** | Frozen Head Pressure : **0-100Kpa**  
Frozen Head Screen : **-12°C to -1°C** | Treatment Handle's Type : **5 Optional (replace applicator)**  
Cooling System : **Semiconductor cooling + water cooling +air cooling**  
Machine Screen : **15.6 inch Collapse screen, 30°-73° adjustable** | Dimension : **41 x 58 x 154cm**



**evolution-med**  
service with a difference



+91-94460-01470  
+91-94460-01472



info@evolutionmed.com  
www.evolutionmed.com